

# PICOSUN™ R-200 ADVANCED

The PICOSUN™ R-200 Advanced ALD systems are suitable for R&D on dozens of applications such as IC components, MEMS devices, displays, LEDs, lasers, and 3D objects such as lenses, optics, jewelry, coins, and medical implants.



## Technical Features

The PICOSUN™ R-200 Advanced ALD system is the global market leader in advanced ALD research tools with hundreds of customer installations. It has become the tool of choice both for companies and research institutes driven by innovation.

The agile design enables the highest quality ALD film depositions together with the ultimate system flexibility to fit future needs and applications. The patented hot-wall design with fully separate inlets and instrumentation enables particle-free processing adaptable on a wide range of materials on wafers, 3D objects, and all nanoscale features. Excellent uniformity even on the most challenging through-porous, ultra-high aspect ratio, and nanoparticle samples is achieved thanks to our proprietary Picoflow™ technology. The PICOSUN™ R-200 Advanced systems are equipped with highly functional and easily exchangeable precursor sources for liquid, gaseous, and solid chemicals. Highly efficient and patented remote plasma option enables deposition of metals without the risk of short-circuiting or plasma damage. Integration with glove boxes, UHV systems, manual and automated loaders, cluster tools, powder chambers, roll-to-roll chambers, and various in situ analytics systems enable efficient and flexible research with good results no matter what your research area is now or might become later on.

### Typical substrate size and type

- 50-200 mm single wafers
- 156 mm x 156 mm solar Si wafers
- 3D objects
- Powders and particles
- Mini-batch
- Porous, through-porous, and high aspect ratio (up to 1:2500)

### Processing temperature

- 50 – 500°C, plasma 450°C (650 °C with heated chuck on request)

### Typical processes

- $\text{Al}_2\text{O}_3$ ,  $\text{TiO}_2$ ,  $\text{SiO}_2$ ,  $\text{Ta}_2\text{O}_5$ ,  $\text{HfO}_2$ ,  $\text{ZnO}$ ,  $\text{ZrO}_2$ ,  $\text{AlN}$ ,  $\text{TiN}$ , metals such as Pt or Ir

### Substrate loading

- Manual loading with a pneumatic lift
- Load lock with magnetic manipulator arm
- Semi-automatic loading with handling robot
- Cassette-to-cassette loading with cluster tools

### Precursors

- Liquid, solid, gas, ozone, plasma
- Up to 12 sources with 6 separate inlets (7 if the plasma option is chosen)


### Options

- Cluster tools, Picoflow™ diffusion enhancer, roll-to-roll chamber, RGA, UHV compatibility,  $\text{N}_2$  generator, gas scrubber, customized designs, glove box integration for inert loading

**Please feel free to contact us for more information or a quotation!**

**The highest quality ALD systems for demanding applications.**


## THE PRINCIPLE OF ALD

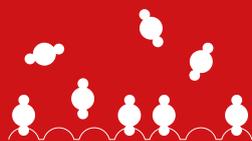

Introduction of molecules  
containing element A.



Adsorption of the molecules  
on the surface.

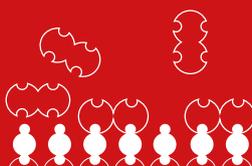

Introduction of molecules  
containing element B and  
reaction with element A on  
the surface.



Completion of one  
monolayer of compound AB.

Repeat cycle till desired film  
thickness is reached.

**picosun**  
AGILE ALD

### Picosun HQ

#### Picosun Oy

Tietotie 3

FI-02150 Espoo, Finland

Tel. +358 50 321 1955

info@picosun.com

#### Picosun Oy

Masalantie 365

FI-02430 Masala, Finland

Tel. +358 50 321 1955

info@picosun.com

### Picosun Branch Offices

#### Picosun Europe GmbH

Tel. +49 1522 449 49 11

sales@picosun.com

#### Picosun USA, LLC

Tel. +1 214 790 0844

Mobile +1 214 490 3951

sales@picosun.com

#### Picosun Asia Pte. Ltd., Singapore

Tel. +65 9756 3265

sales@picosun.com

#### Picosun Taiwan Co. Ltd.

Tel. +886 90 515 2985

sales@picosun.com

#### Picosun China Co. Ltd.

Tel. +358 40 480 3449

sales@picosun.com

#### Picosun Japan Co. Ltd.

Tel. +81 3 6431 9500

Mobile +81 90 5198 8131

sales@picosun.com