

# PICOSUN™ R-200 STANDARD

The PICOSUN™ R-200 Standard ALD systems are suitable for R&D on dozens of applications such as IC components, MEMS devices, displays, LEDs, lasers, and 3D objects such as lenses, optics, jewelry, coins, and medical implants.



## Technical Features

The PICOSUN™ R-200 Standard ALD system is the market leader in thermal ALD research tools. It has become the tool of choice both for companies and research institutes driven by innovation.

The agile design enables the highest quality ALD film depositions together with the ultimate flexibility of the system to fit future needs and applications. The patented hot-wall design with fully separate inlets and instrumentation enables particle-free processing adaptable on a wide range of materials on wafers, 3D objects, and all nanoscale features. Excellent uniformity even on the most challenging through-porous, ultra-high aspect ratio, and nanoparticle samples can be achieved thanks to our proprietary Picoflow™ technology. The PICOSUN™ R-200 Standard systems are equipped with highly functional and easily exchangeable precursor sources for liquid, gaseous, and solid chemicals. Integration with glove boxes, powder chambers, and various in situ analytics systems enable efficient and flexible research with good results, no matter what your research area is now, or what it might become later on.

### Typical substrate size and type

- 50-200 mm single wafers
- 156 mm x 156 mm solar Si wafers
- 3D objects
- Powders and particles
- Mini-batch
- Porous, through-porous, and high aspect ratio (up to 1:2500)

### Processing temperature

- 50 – 500°C

### Typical processes

- $\text{Al}_2\text{O}_3$ ,  $\text{TiO}_2$ ,  $\text{SiO}_2$ ,  $\text{Ta}_2\text{O}_5$ ,  $\text{HfO}_2$ ,  $\text{ZnO}$ ,  $\text{ZrO}_2$ ,  $\text{TiN}$ ,  $\text{AlN}$ , metals such as Pt or Ir

### Substrate loading

- Manual loading with a pneumatic lift
- Load lock with magnetic manipulator arm

### Precursors

- Liquid, solid, gas, ozone
- Up to 6 sources with 4 separate inlets


### Options

- Picoflow™ diffusion enhancer, RGA,  $\text{N}_2$  generator, gas scrubber, customized designs, glove box compatibility for inert loading

The highest  
quality ALD  
systems for  
demanding  
applications.

Please feel free to contact us for more  
information or a quotation!


## THE PRINCIPLE OF ALD

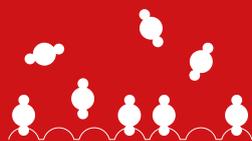

Introduction of molecules  
containing element A.



Adsorption of the molecules  
on the surface.

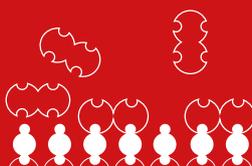

Introduction of molecules  
containing element B and  
reaction with element A on  
the surface.



Completion of one  
monolayer of compound AB.

Repeat cycle till desired film  
thickness is reached.

**picosun**  
AGILE ALD

### Picosun HQ

#### Picosun Oy

Tietotie 3  
FI-02150 Espoo, Finland  
Tel. +358 50 321 1955  
info@picosun.com

#### Picosun Oy

Masalantie 365  
FI-02430 Masala, Finland  
Tel. +358 50 321 1955  
info@picosun.com

### Picosun Branch Offices

#### Picosun Europe GmbH

Tel. +49 1522 449 49 11  
sales@picosun.com

#### Picosun USA, LLC

Tel. +1 214 790 0844  
Mobile +1 214 490 3951  
sales@picosun.com

#### Picosun Asia Pte. Ltd., Singapore

Tel. +65 9756 3265  
sales@picosun.com

#### Picosun Taiwan Co. Ltd.

Tel. +886 90 515 2985  
sales@picosun.com

#### Picosun China Co. Ltd.

Tel. +358 40 480 3449  
sales@picosun.com

#### Picosun Japan Co. Ltd.

Tel. +81 3 6431 9500  
Mobile +81 90 5198 8131  
sales@picosun.com